

Dental Instruments



AN

Endodontics

Endodoncia

Endodontie

Quality makes Dentistry Better

 **ACUMED**
International Co.

Serving Healthcare Worldwide Since 1984

Classic

*Satin Finish
8mm Hollow Handle
Fully Autoclavable*

**ULTRA-LIGHT
Stainless Steel Made**


DID YOU KNOW?

ACUMED Endodontic Explorers has, Long tip with sharp point. They are used to locate the opening of small canals during endodontic procedures. All Acumed endodontic explorers are made of spring tempered stainless steel. Due to their flexibility they can retain their shape and resist breakage.

All products are available with variety of handle designs, if you required some specific handle for your products, dont hesitate to ask.

Quality makes Dentistry Better

All ACUMED hand instruments are available in given handle designs. Just add the handle code after article number to order special handles.

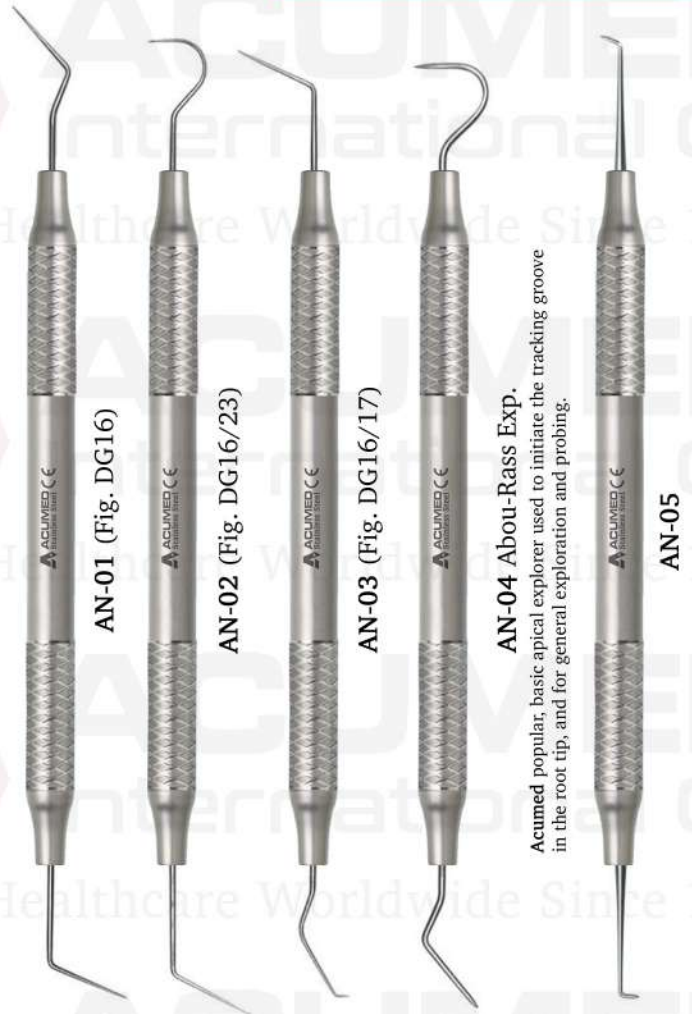

ACUMED Endodontic Excavators have extra-long terminal shanks to reach deep into the cavity preparation. All are made of high quality non rusty stainless steel that is formed and precision ground by expert craftsmen, then hardened for the ultimate in sharp edge retention and durability.

Classic

*Satin Finish
8mm Hollow Handle
Fully Autoclavable*



**ULTRA-LIGHT
Stainless Steel Made**




Endodontic

Acumed Endodontic Plugger tips are made from spring tempered stainless steel wire for maximum durability and used to vertically condense gutta-percha within the canal.

21mm Endodontic Luks

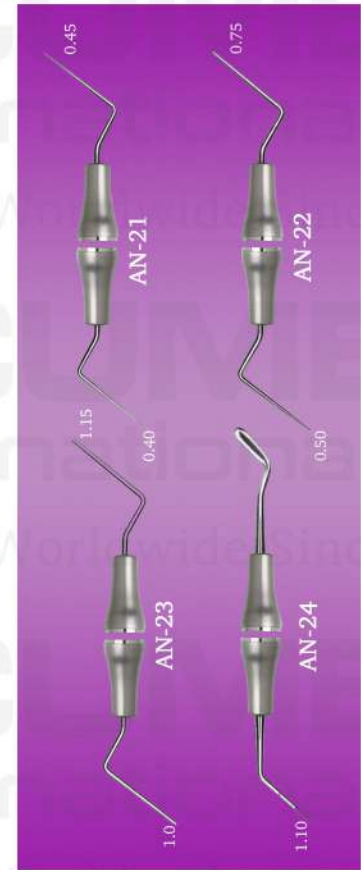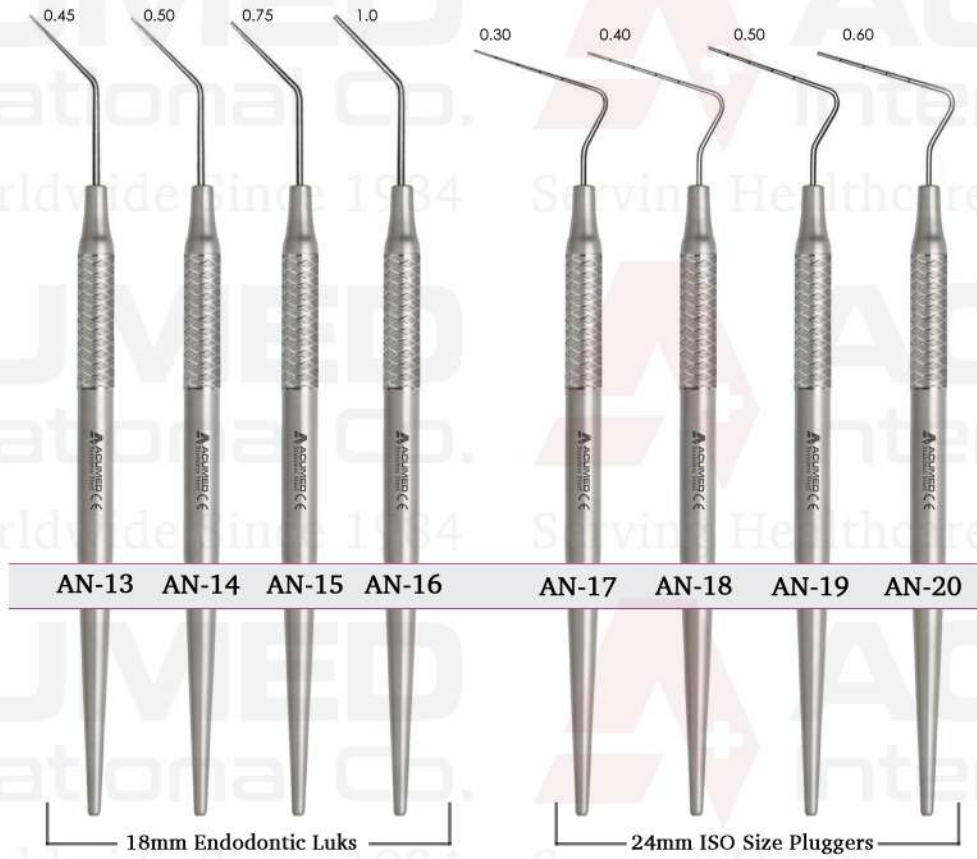

Used to compact filling material during vertical condensation

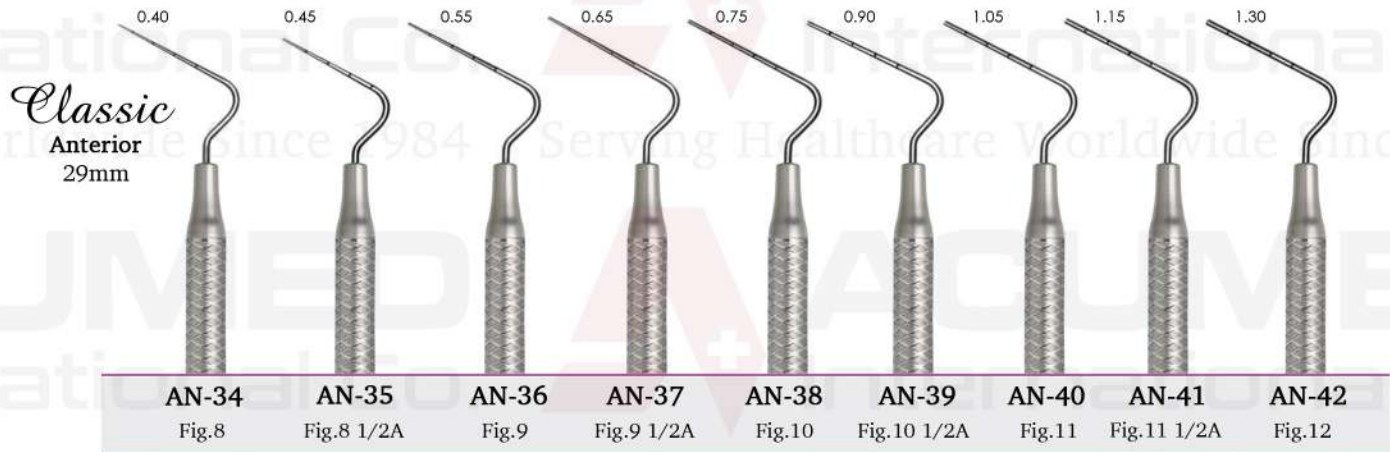

ACUMED hand instruments with a tapered, pointed end, used to condense and compress materials such as gutta percha in root canals.

See Instructions before use.

sales@acumedinternational.com.

ISO Sized Spreaders

AN-43 (24mm)
AN-44 (24mm)
AN-45 (28mm)
AN-46 (28mm)

ACUMED
Serving Healthcare Worldwide Since 1984
Quality makes Dentistry Better

ISO9001, ISO13485, CE Certified Company
Jinnah Town, Capital Road, Siakot - 51310, Pakistan.

Heat Carriers



AN-47 Fig.3 - 28mm AN-48 Fig.GP2 - 25mm AN-49 Fig.0 - 25mm AN-50 Fig.00 - 25mm AN-51 Fig.MA57 - 27mm AN-52 Fig.GP3 - 30mm




AN-53 GP1 - 20mm AN-54 1S Wakal - 22mm AN-55 2S Wakal - 22mm AN-56 D11 Stand - 23mm AN-57 D11S Thin - 23mm AN-58 D11T Thinner - 22mm AN-59 D11TS Thinner - 22mm

TRUE BLACK TITANIUM

Used to compact filling material during lateral condensation. These instruments are made of nickel titanium.

ULTRA-LIGHT Stainless Steel Made

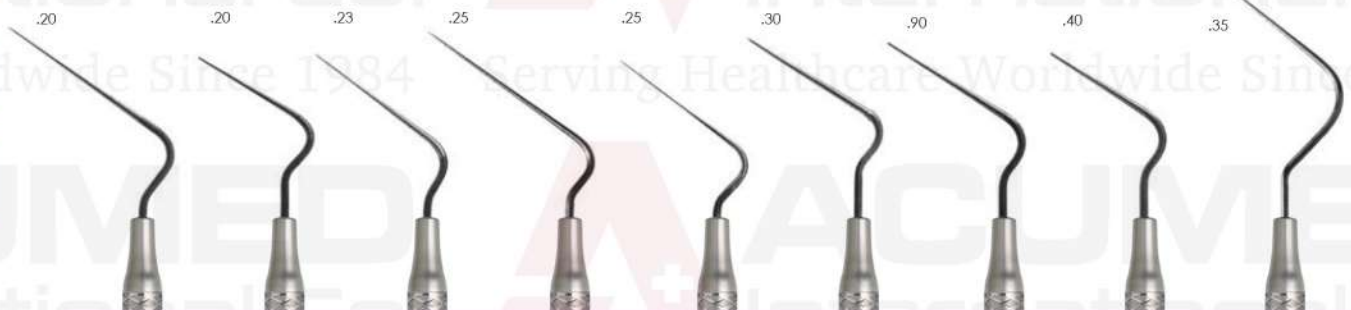

AN-60 Fig. MA57 AN-61 Fig. D11TS AN-62 Fig. D11T AN-63 Fig. 25 AN-64 Fig. 255 AN-65 Fig. D11 AN-66 Fig. 40 AN-67 Fig. 40S AN-68 Fig. 3NT

Satin Finish
8mm Hollow Handle
Fully Autoclavable

Classic

Perforation repair instruments and Appical pluggers for compact filling material during vertical condensation.



Micro Endodontic Burnishers / Pluggers
Used to condense, smooth and polish amalgam



Satin Finish
8mm Hollow Handle
Fully Autoclavable



Retro Felling Instruments - used to work with ultrasonic preparation tips



Classic

Ruddle Currettes are used to make a sharp dissection of the periosteal tissue, thus preventing tissue damage. These relatively small currettes allow for undermining elevation of the flap, which along with the micro-surgical incision are the first two steps in automatic flap management.



ACUMED Endodontic Currettes are extremely sharp and able to cut through attaching fibers with the efficiency of a scalpel. This allows the skilled clinician to provide patients with the shortest recovery time possible. All of these small currettes are manufacturered from highest-grade surgical stainless steel and are autoclavable.

ACUMED Endo currettes are designed for curettage, cyst removal and tooth socket debridement.



Pluggers and Burnishers

Pluggers are used to condense filling materials into cavity preparation. Ends are serrated. Burnishers are used to condense, smooth and polish amalgam



*Satin Finish
8mm Hollow Handle
Fully Autoclavable*



Classic

Endodontic Instrument Retrieving Forceps
Remove broken endo instruments
Tips are 20mm long by 0.8mm diameter overall, with 5mm serrated edges at the very tip of the instrument.



AN-98
Diamond Coated Tips,
Gold Plated Rings

Diamond Peets Forceps are excellent for removing separated files or posts. The thin tips are fine enough to access tight occlusal openings, while the diamond coated beaks provides unsurpassed gripping performance.




PEET'S Forceps
CURVED HANDLES,
30° Angled,
12cm

BEST SELLER



AN-100
Diamond Coated Tips,
Gold Plated Rings

Diamond Peets Forceps are excellent for removing separated files or posts. The thin tips are fine enough to access tight occlusal openings, while the diamond coated beaks provides unsurpassed gripping performance.

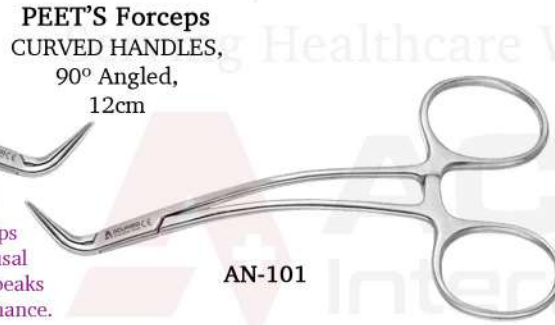

PEET'S Forceps
CURVED HANDLES,
90° Angled,
12cm

AN-101



Nerve canal pliers 15cm



AN-102 **AN-103**

Endodontic Tweezers



AN-104 No LOCK **AN-105** with LOCK **AN-106** No LOCK



AN-107 Stieglitz Root Forceps, Tapered Tips, Straight, 14cm **AN-108** Stieglitz Root Forceps, Tapered Tips, 45°, 13.5cm **AN-109** Stieglitz Root Forceps, Tapered Tips, 90°, 12cm

BEST SELLER

For root fragments For upper root and teeth For upper root and teeth



AN-110 AN-111 AN-112 AN-113



Straight
AN-128

Upward
AN-129

Left
AN-130

Right
AN-131

The complete bendable MTA carrier syringe and plunger for applying MTA materials. The three piece design allows for replacing sections. The syringe, plunger, and body of the MTA carrier can be purchased separately.



- AN-114** MTA Carrier 0.8mm Bendable
- AN-115** MTA Carrier 0.99mm Bendable
- AN-116** MTA Carrier 1.6mm Bendable

Endodontic Aspirators



- AN-117** 0.5mm Micro Suction Tip
- AN-118** 1.0mm Micro Suction Tip
- AN-119** 2.0mm Micro Suction Tip

- AN-120** 1.0mm Micro Suction Tip w/Tail
- AN-121** 2.0mm Micro Suction Tip w/Tail

- AN-122** 1.7mm Long Endo Aspirator
- AN-123** 2.3mm Long Endo Aspirator

- AN-124** 1.7mm Long Endo Aspirator
- AN-125** 2.3mm Long Endo Aspirator



- AN-126** 1.0mm Cleaning Stylet
- AN-127** 2.0mm Cleaning Stylet

Super Endo Dental Buchanan Hand Pluggers

Condensers are available individually or in a full set of 4 condensers.

One end is made of nickel titanium, can be used for treatment in curved canals.

The other end, flat and made of stainless steel, is ideal for compaction in the coronal and middle thirds of canal.

Handles are colour coded for quick identification.

25mm Long Tips



Specification :

Code: AN-SDB-1 #0Stainless tip size.02/0.75 NiTi Tip Size.03/0.25 Yellow Tip

Code: AN-SDB-2 #1Stainless tip size.02/0.8NiTi Tip Size.03/0.4 Red Tip

Code: AN-SDB-3 #2Stainless tip size.02/.1.3NiTi Tip Size.03/0.7 Blue Tip

Code: AN-SDB-4 #3Stainless tip size.02/.68NiTi Tip Size.03/0.09 White Tip

EXTREMELY DELICATE - HARD - POWERFUL GRIP

The diamond-like coating enables a safe grip for smallest parts as e. g. endodontic instruments, tiny root fragments, stuck matrix strips etc.

Manufactured of high quality, medical stainless steel.



AN-ENDO1

AN-ENDO2

Carr Endodontic Retractors

Used to hold mucoperiosteal flaps.

Retracts cheeks, lips, and tongue.

Serrated tips for reduce Slippage.

Available in 45 Degree and 90 Degree angles.

AN-CARR45

AN-CARR90





Acumed Logo and Name are Registered and Trade Mark of Acumed



Scan our mobile website to your smartphone for quick reference

Acumed International Co.
Corporate Head
Jinnah Town, Capital Road,
Sialkot - 51310, Pakistan.

Email: sales@acumedco.com

Website: www.acumedco.com